

# Dr. Yoshi Tips




*Yoshi Terauchi*

- Received DDS in 1993
- Completed residency at Tokyo Medical & Dental University in 1995, where he also received PhD from the Department of Endodontics and internationally
- Published several articles in peer-reviewed journals nationally and internationally
- Authored in several chapters in textbooks including the 11th and the 17th edition of "Pathways of the Pulp" and "Endodontics: the 6th edition of Principles and Practice"
- Lectured nationally and internationally and was exposed twice on National TV for modern endodontics.

## E87/ED87

(Jointly Developed with Dr. Yoshi Terauchi)

Tip for broken file removal.  $\Phi 0.1$ mm tip head, extremely thin and conical, used to remove dentin around the broken file.

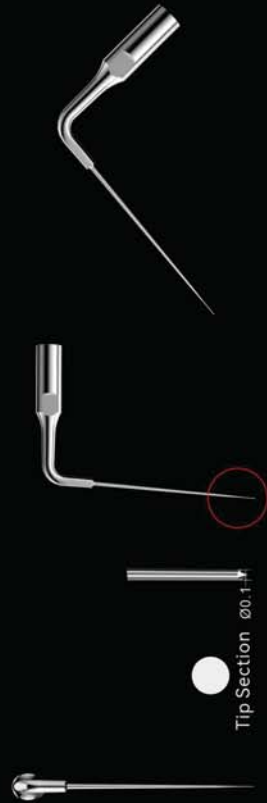


Tip Section  $\Phi 0.1$

## E88/ED88

(Jointly Developed with Dr. Yoshi Terauchi)

Tip for broken file removal. 0.1mm thickness tip head, extremely-thin sword shaped, front and back directed, used to remove dentin around the broken file.

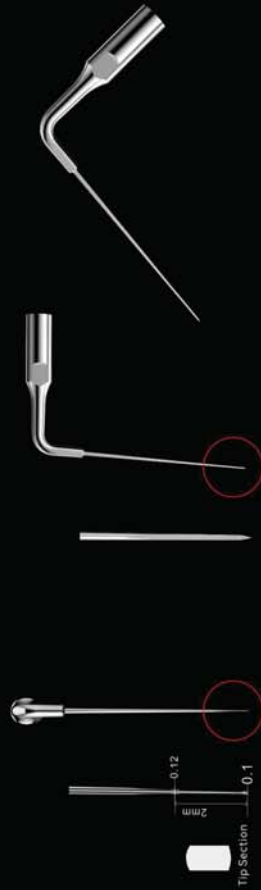


Tip Section  
2mm  
0.1

## E89/ED89

(Jointly Developed with Dr. Yoshi Terauchi)

Tip for broken file removal. 0.1mm thickness tip head, extremely-thin sword shaped, left and right directed, used to remove dentin around the broken file.

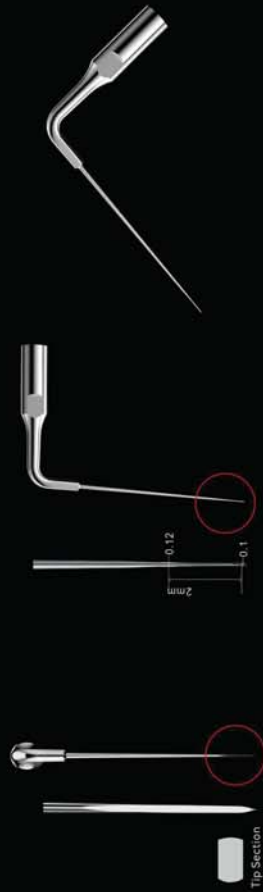

Tip Section  
2mm  
0.1